

Eisner™

A division of Multigate
Medical Products Pty Ltd



Insulated Tubes, Sheaths & Locking Nuts

Eisner™ laparoscopic
instruments deliver
unmatched durability,
comfort and efficiency
to support a seamless
performance every time.

1993

POWERING
LAPAROSCOPIC
SURGERY
SINCE 1993



DESIGN

Our instruments are designed for optimal balance, feel, comfort and responsiveness to deliver a seamless connection between the hand and the surgical site.

SAFETY & RELIABILITY

The craftsmanship and the quality of materials we select combine to deliver high performance, electro-safe instruments that deliver unrivalled durability, longevity and reliability.



• **Easy to re-process** Our instruments are designed to simplify re-processing.

• **Clinical performance** Our hand-held instruments are modular allowing you to build cost effective instrument trays whilst maximising patient safety and clinical performance.




INSULATED TUBES, SHEATHS & LOCKING NUTS

Our insulated tubes are designed for easy re-processing and to reduce the potential for fluids, debris and bio-burdens to enter under the insulation.



10MM INSULATED TUBE

Durable Insulation material density increased three-fold to minimise splitting, cracking or deformity.

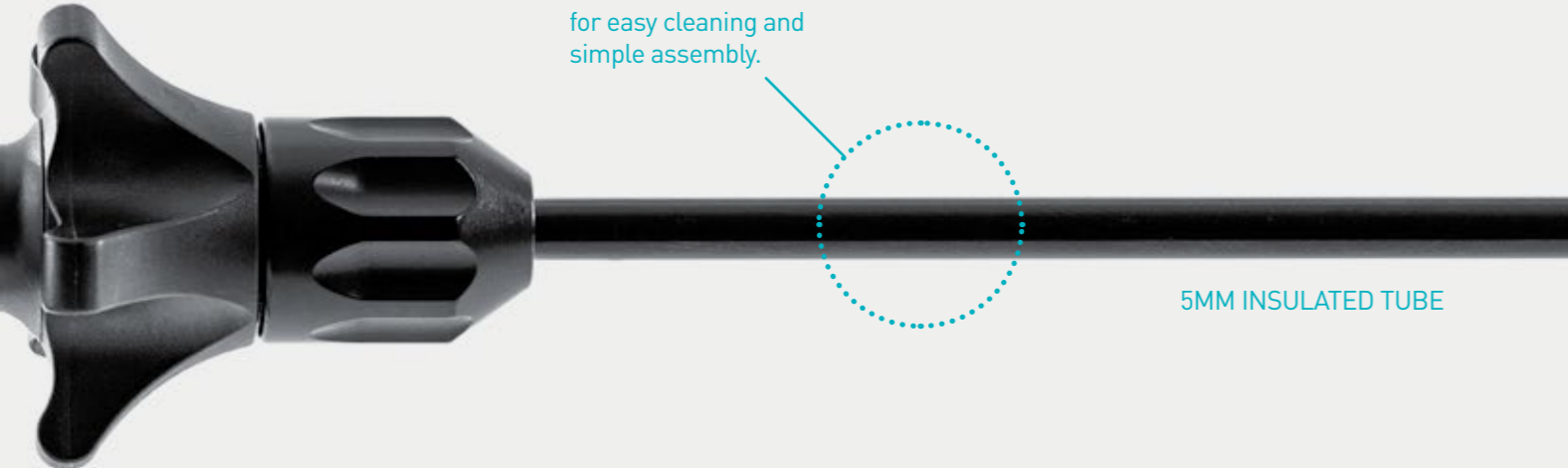

5MM INSULATED TUBE

Ease Our Ezi-Clean™ Insulated Tubes have a removable sheath for easy cleaning and simple assembly.




Our range of coloured locking nuts help you to code your instruments trays to specific laparoscopic procedures.


Insulated Tubes, Sheaths & Locking Nuts (5mm)






Our insulated tubes are made from surgical grade stainless steel and have a high performance PPSU polymer insulation designed to resist heat, impact and chemicals.

Tubes and sheaths are available in standard and bariatric lengths and in both 5mm and 10mm diameters.

Safety-Lip	DESCRIPTION	RE-ORDER 33cm	RE-ORDER 45cm	FEATURES	MIN QTY
	5mm Safety-Lip Insulated Tube with locking nut	IT-05/33	IT-05/45	Steel Safety-lip Sealed Tip	1






Ezi-Clean™	DESCRIPTION	RE-ORDER 33cm	RE-ORDER 45cm	FEATURES	MIN QTY
	5mm Ezi-Clean™ Insulated Tube with locking nut	IT-05/33.EZC	IT-05/45.EZC	Insulation Tube is Removable Re-usable & Replaceable	1

	DESCRIPTION	RE-ORDER	RE-ORDER	FEATURES	MIN QTY
	Black Insulation Sheath Only (Pack of 5)	EZC-33	EZC-45	Replaceable Insulation Sheaths	1 Pack

Locking Nuts 5mm	DESCRIPTION	RE-ORDER	SIZE	FEATURES	MIN QTY
	5mm Eisner Locking Nut - Standard Black	ELN-05SB	Ø5mm	Standard Black	1
	5mm Eisner Locking Nut - Blue	ELN-05B	Ø5mm	Blue	1
	5mm Eisner Locking Nut - Green	ELN-05G	Ø5mm	Green	1
	5mm Eisner Locking Nut - Red	ELN-05R	Ø5mm	Red	1
	5mm Eisner Locking Nut - Yellow	ELN-05Y	Ø5mm	Yellow	1

Insulated Tubes, Sheaths & Locking Nuts (10mm)

Molded Hardcoat	DESCRIPTION	RE-ORDER 33cm	RE-ORDER 45cm	FEATURES	MIN QTY
	10mm Molded Hardcoat Insulation Tube with Locking Nut	IT-10/33	IT-10/45	Thickened Hardcoat Insulation	1

Locking Nuts 10mm	DESCRIPTION	RE-ORDER	SIZE	FEATURES	MIN QTY
	10mm Eisner Locking Nut - Standard Black	ELN-10SB	Ø10mm	Standard Black	1
	10mm Eisner Locking Nut - Blue	ELN-10B	Ø10mm	Blue	1
	10mm Eisner Locking Nut - Green	ELN-10G	Ø10mm	Green	1
	10mm Eisner Locking Nut - Red	ELN-10R	Ø10mm	Red	1
	10mm Eisner Locking Nut - Yellow	ELN-10Y	Ø10mm	Yellow	1

We work at the intersection between the surgeon and their patient, designing high performance, responsive laparoscopic instruments that offer a seamless connection between the hand and the surgical site.



SUPPORT

With Eisner™ you don't get a call centre, just a dedicated team of specialists ready to tailor your Eisner™ experience.

Experience the difference

Book an obligation free trial

Our team can organise an obligation free, trial assessment of your instruments so you can see and feel the difference before you buy.

Trade-in / replacement program

Manage your tray upgrades

If you are considering decommissioning an instrument our trade-in/replacement service will assist you to address one off replacements or the gradual upgrade of instrument sets. We accept Eisner™ and non-Eisner™ instruments.

Care

Optimise Your Tray Performance

When you choose Eisner™ products our team are available to provide free comprehensive in-service training, support and tray assessments to guide servicing and repair.





Enhance your every move

Warranty | Eisner™ offers a 12 month warranty against manufacturing defect on all metal components of the equipment. Damage and / or breakages caused by improper use or abuse of the equipment are not covered by warranty. Normal wear and tear is not covered. Warranty is null and void where components are mixed with non-Eisner™ instruments.

Service & Repairs | Instruments returned to us for repair/replacement will be assessed and a written quotation provided for you to make a decision on further action. In most cases for instruments beyond repair we can offer a trade-in replacement. Eisner™ keeps factory original spare parts in stock. All repairs should be returned to Eisner™ for repair in order to preserve the integrity of your instruments.

Recommended Lifespan of Instruments | Whilst Eisner™ laparoscopic instruments will function for many years, we recommend replacement prior to any obvious signs of wear and stress which could lead to an issue during surgery. We recommend a lifespan of 3 years on our handles and inserts and suggest that after that time, they are traded in on a new component. Insulated Tubes: By virtue of the fact that insulated tubes will be replaced once the insulation is damaged or compromised in any way, it is not envisioned that the lifespan of a sheath will be beyond 12 months. Instruments should be taken out of circulation immediately if it is found that identifying labelling is no longer visible.

EisnerTM

A division of Multigate
Medical Products Pty Ltd

Pacifica Business Park
Unit 18, 97 Canterbury Road
Kilsyth VIC 3137
Australia
T +61 (0)3 9728 4770
F +61 (0)3 9728 6770
eisner@multigate.com.au
www.eisner.com.au

For an obligation
free trial or for more
information call
us today on:
(03) 9728 4770
eisner.com.au